



**WHYNTER Bamboo Bathroom Scale**

# **OWNER'S MANUAL**

**MODEL: BH-2000**

# WHYNTER Bamboo Electronic Bathroom Scale

Your new electronic personal scale has been designed to accurately indicate your weight gain or loss over a period of time and should give many years of service under normal use.

## SPECIFICATIONS

- Equipped with a high precision strain gauge sensor system
- Capacity 150kgs(330 lbs/24 st)
- Division 100g/0.2lb
- Foot-tap switch
- Automatic zero resetting
- Automatic switch off
- Low battery / Over-load indicator

## POWER SUPPLY

- battery : □2 x CR2016 □1 x 9 V □1 X CR2032 □2 X CR2032
- For lithium battery powered scale, please pull the insulation tape out before first using.



- When the LCD indicate “Lo”, open the battery door and replace

the battery.



## OPERATION


STEP 1: Put the scale on a hard & flat surface (avoid carpet or soft surface).

STEP 2: Switch on the scale by slightly tapping on the platform, and then wait a few seconds, the following will appear.



(NOTE: Please make sure the "0.0" appears on the LCD display before standing on the scale.)

STEP 3: Gently step on the platform, put your feet evenly on the scale without moving. The scale will compute your weight and it will be displayed in LCD display.

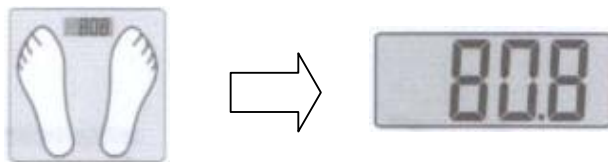

## AUTOMATIC SWITCH OFF

When you step off the scale, it will automatically switch off in approximately 6 seconds.

## OVER LOAD INDICATOR

When the scale is over loaded (over 150kg/ 330lbs/ 24st), "Err" will

appear on the screen.



### ADVICE FOR USE AND CARE

- Always place your scale on a hard & flat surface to ensure the greatest weighing accuracy.
- Weigh yourself without clothing and footwear, before meals and at the same time of day. Morning is the best time to weight yourself.
- Your scale rounds up or down to the nearest 100g, if you weigh yourself twice and get tow different readings, your weight lies between the two readings.
- Clean your scale with a damp cloth but do not allow water to enter the scale. Do not use chemical cleaning agents.
- Do submerge the scale in water.
- Do not store the scale in an upright position.
- Remove the batteries if the scale is not being used for a long time.
- Do not drop the scale or jump on it.
- Do not use for commercial trade.

## **180 DAYS LIMITED WARRANTY**

This WHYNTER scale is warranted, to the original owner within the 48 continental states, for 180 days from the date of purchase against defects in material and workmanship under normal use and service. Should your WHYNTER scale prove defective within 180 days from the date of purchase, return the defective part or unit, freight prepaid (within one month of purchase; after one month to one year, customer will be responsible for freight cost to Whynter's service department), along with an explanation of the problem/issue. Please carefully package the scale in its original packaging material to avoid damage in transit. Please retain original box and packaging material. Under this warranty, Whynter will repair or replace any parts found defective. This warranty is not transferable. After the expiration of the warranty, the cost of labor and parts will be the responsibility of the original owner of the unit.

### THIS WARRANTY DOES NOT COVER:

- Acts of God, such as fire, flood, hurricanes, earthquakes and tornadoes.
- Damage, accidental or other wise, to the scale while in the possession of a consumer not caused by a defect in material or workmanship;
- Damage caused by consumer misuse, tampering, or failure to follow the care and special handling provisions in the instructions.
- Damage to the finish of the case, or other appearance parts caused by wear.
  - 1) Surface
  - 2) Damage caused by repairs or alterations of the scale by anyone other than authorized by the manufacturer.
  - 3) Freight and Insurance cost for the warranty service.

To obtain service or information, contact

Whynter Innovations Group via Email at [service@whynter.com](mailto:service@whynter.com) or  
call 866-WHYNTER.